**3M Unitek
Orthodontic Products**

2724 South Peck Road
Monrovia, CA 91016 USA

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Unitek™ Temporary Anchorage Device (TAD)
Constant Force Coil Springs

Instructions for Use

ENGLISH

Indications for use

This product is intended for use in orthodontic treatment.


WARNINGS:

Single use only. Do not reuse. Reuse may result in compromising product integrity or may create risk of patient infection or cross contamination.

This product contains nickel and/or chromium. A small percentage of the population is known to be allergic to nickel and/or chromium. If an allergic reaction occurs, direct the patient to consult a physician.

Locking Eyelet

Hook Eyelet

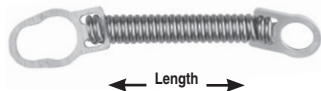Product description

There are three components to the preassembled Unitek™ Temporary Anchorage Device (TAD) Constant Force Coil Springs: the locking eyelet, the hook eyelet, and the coil spring. The Unitek TAD Constant Force Coil Springs are preassembled at both ends with stainless steel eyelets that allow for easy engagement of various appliance hooks, as well as the Unitek™ Temporary Anchorage Device (TAD). The eyelet that engages the TAD has a specially designed locking system that helps keep it attached under the most severe conditions. The coil springs are made of nickel-titanium. The locking eyelet is a constant length of 4.7 mm. The hook eyelet is a constant length of 3.3 mm. The spring working length varies by size.

Part Number	Length	Force	Target Load (g)	Maximum Activation Length
504-2221	1.75 mm	Light	150	8 mm
504-2222		Medium	200	8 mm
504-2223		Heavy	250	8 mm
504-2231	3 mm	Light	150	12 mm
504-2232		Medium	200	12 mm
504-2233		Heavy	250	12 mm
504-2261	6 mm	Light	150	18 mm
504-2262		Medium	200	18 mm
504-2263		Heavy	250	18 mm
504-2281	8 mm	Light	150	22 mm
504-2282		Medium	200	22 mm
504-2283		Heavy	250	22 mm

Placement Instructions

First place the Hook Eyelet over the hook attached to the orthodontic appliances. Attaching the bracket hook first minimizes excessive force being placed on the TAD. Then, gripping the leading edge of the smaller diameter hole of the Locking Eyelet with Weingart pliers, pull gently to slip the larger diameter hole of the Locking Eyelet over the TAD and release the eyelet. The smaller hole of the eyelet will then positively engage the neck of the TAD and lock the spring in place.

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3M and Unitek are trademarks of 3M or
© 2014 3M. All rights reserved. REF 011-641-2ML 1410